

SUBSTITUTE

SPECIFICATION

METHOD AND APPARATUS FOR RECONSTRUCTING A LIGAMENT

Reference To Pending Prior Patent Application

CROSS-REFERENCE TO RELATED APPLICATIONS

[0001] This patent application claims benefit of pending prior U.S. Provisional Patent Application Serial No. 60/237,817, filed 10/03/00 by E. Marlowe Goble for METHOD AND APPARATUS FOR RECONSTRUCTING A LIGAMENT, which patent application is hereby incorporated herein by reference.

BACKGROUND OF THE INVENTION

1. The Field Of The of the Invention

[0002] This invention relates to medical procedures and apparatus in general, and more particularly to bone fixation systems and related medical procedures and apparatus for that can be used in reconstructing a ligament.

Background Of The Invention2. The Relevant Technology

[0003]—A ligament is a piece of fibrous tissue which connects one bone to another.

[0003] Ligaments are frequently damaged (e.g., detached or torn or ruptured, etc.) as the result of injury and/or accident. A damaged ligament can impede proper motion of a joint and cause significant pain.

[0004] Various procedures have been developed to repair or replace a damaged ligament. The specific procedures used depend on the particular ligament which is to be restored and on the extent of the damage.

Page 1

[0005] One ligament which is frequently damaged as the result of injury and/or accident is the anterior cruciate ligament (ACL). Looking now at Fig. 1, the ACL 5 extends between the top of the tibia 10 and the bottom of the femur 15. A damaged ACL can cause instability of the knee joint and cause substantial pain and arthritis.

[0006] Numerous procedures have been developed to restore the ACL through a graft ligament replacement. In general, and looking now at Fig. 2, these ACL replacement procedures involve drilling a bone tunnel 20 through tibia 10 and up into femur 15. Then a graft ligament 25, consisting of a harvested or artificial ligament or tendon(s), is passed through the tibial portion 30 of tunnel 20 (sometimes referred to as the "tibial tunnel"), across the interior of the joint, and up into the femoral portion 35 of tunnel 20 (sometimes referred to as the "femoral tunnel"). Then a distal portion of graft ligament 25 is secured in femoral tunnel 35, and a proximal portion of graft ligament 25 is secured in tibial tunnel 30.

BRIEF DESCRIPTION OF THE DRAWINGS

[0007] These and other objects and features of the present invention will be more fully disclosed or rendered obvious by the following detailed description of the preferred embodiments of the invention, which are to be considered together with the accompanying drawings wherein like numbers refer to like parts and further wherein:

Fig. 1 is a schematic side view of a knee joint, showing an ACL extending between the top of the tibia and the bottom of the femur;

Fig. 2 is a schematic side view of the same knee joint, showing portions of an ACL reconstruction;

Fig. 3 is a schematic front view of a knee joint, illustrating a drill guide inserted into a bone tunnel;

Figs. 4-6 are schematic views of one preferred form of a drill guide;

Figs. 7-12 are schematic views illustrating an ACL reconstruction procedure;

Figs. 13 and 14 are schematic views illustrating an ACL reconstruction procedure with bone blocks attached to a tendon; Figs. 15 and 16 are schematic views illustrating an alternative ACL reconstruction procedure with a wire suture having a loop at its distal end; and

Figs. 17-21 are schematic views of another alternative reconstruction procedure with a cannulated screw and associated threaded collet.

[0008] Figure 1 is an elevated side view of a knee joint having an anterior cruciate ligament;

[0009] Figure 2 is a side view of the knee joint shown in Figure 1 with the anterior cruciate ligament replaced by a graft ligament;

[0010] Figure 3 is a cross sectional side view of a knee joint having a drill guide partially disposed within a bone tunnel formed on the tibia and femur;

Page 3

- [0011] Figure 4 is a perspective view of an alternative embodiment of the drill guide shown in Figure 3;
- [0012] Figure 5 is a top perspective view of a portion of the drill guide shown in Figure 4;
- [0013] Figure 6 is a bottom perspective view of a portion of the drill guide shown in Figure 4;
- [0014] Figure 7 is a cross sectional side view of the knee joint shown in Figure 3 having the drill and barrel guide removed;
- [0015] Figure 8 is a cross sectional side view of the knee joint shown in Figure 7 with a suture engaged with the endosteal guide;
- [0016] Figure 9 is a cross sectional side view of the knee joint shown in Figure 8 with a cannulated screw passed over the suture and driven into the angled tunnel;
- [0017] Figure 10 is a cross sectional side view of the knee joint shown in Figure 9 with the endosteal guide removed from the bone tunnel;
- [0018] Figure 11 is a cross sectional side view of the knee joint shown in Figure 10 with a graft ligament attached to the suture;
- [0019] Figure 12 is a cross sectional side view of the knee joint shown in Figure 11 with the draft ligament drawn into the bone tunnel and a crimp secured to the suture;
- [0020] Figure 13 is a cross sectional side view of the knee joint shown in Figure 12 wherein the graft ligament is depicted as a bone-tendon-bone graft that is coupled with the suture and is disposed outside of the bone tunnel;
- [0021] Figure 14 is a cross sectional side view of the knee joint shown in Figure 13 wherein the bone-tendon-bone graft is drawn into the bone tunnel;
- [0022] Figure 15 is a cross sectional side view of the knee joint shown in Figure 7 wherein a grasper disposed within the bone tunnel is grasping a suture passed through the angled tunnel;

- [0023] Figure 16 is a cross sectional side view of the knee joint shown in Figure 7 wherein the suture has been pulled through the bone tunnel by the grasper;
- [0024] Figure 17 is an elevated side view of an assembly including a suture, cannulated screw, collet, and drivers;
- [0025] Figure 18 is a back perspective view of the cannulated screw and collet shown in Figure 17;
- [0026] Figure 19 is a cross sectional side view of the cannulated screw shown in Figure 18;
- [0027] Figure 20 is a cross sectional side view of the assembled suture, cannulated screw, collet shown in Figure 17; and
- [0028] Figure 21 is a perspective view of the driver for driving the cannulated screw and the driver for driving the collet as shown in Figure 17.

DETAILED DESCRIPTION OF THE PREFERRED EMBODIMENTS

[0038][0029] There are numerous ways in which graft ligament 25 may be loaded into bone tunnel 20 and then secured in position. The present invention is directed to a new method for positioning a graft ligament 25 in bone tunnel 20 and for securing the graft ligament in position, and to new apparatus for use in the same.

[0039][0030] More particularly, and looking now at Fig. 3, after bone tunnel 20 (consisting of tibial tunnel 30 and femoral tunnel 35) have has been drilled, a drill guide 40 is inserted into the bone tunnel, i.e. 20. Specifically, the drill guide 40 has its an endosteal guide 45 passed up tibial tunnel 30, across the interior of the knee joint, and then up femoral tunnel 35. Then an anterolateral portal is made on the anterior border of the IT band, and an angled tunnel 50 is drilled from the periosteum to the proximal end of femoral tunnel 35. Thus, angled tunnel 50 opens on, and communicates with, femoral tunnel 35. In this respect it should be appreciated that drill guide 40 is essentially a cross-pin drill guide of the sort well known in the art, except modified so as to enable the angled tunnel 50 to be drilled at an acute angle (e.g., 45 degrees) to the axis of femoral tunnel 35.

[0040][0031] If desired, drill guide 40 can be set so as to always drill the angled tunnel 50 at the same angle (e.g., 30 or 45 or 60 degrees, etc.), or the approach may be straight, i.e., aligned with the femoral tunnel (e.g., at 0 degrees), or drill guide 40 can be configured so as to be adjustable, whereby the surgeon can select the specific angle of angled tunnel 50. By way of example but not limitation, some details regarding the construction of one preferred form of drill guide 40 is shown in Figs. 4-6.

<u>[10041][10032]</u> Looking next at Fig. 7, after angled tunnel 50 is drilled, the drill and barrel guide are removed from drill guide 40, and the drill guide is . The drill guide 40 is then rotated about the axis of bone tunnel 20 so as to swing the drill guide's head of drill guide 40 away from the

Page 6

entrance to angled tunnel 50, whereby to provide convenient access to the entrance of angled

tunnel 50.

[0042][0033] Looking next at Fig. 8, a braided wire suture 55 is inserted down angled tunnel 50

and "grabbed" by the distal end of endosteal guide 45. To this end, the endosteal guide 45 may

include some sort of grippers at its distal end, or the end of the wire suture 55 may include a loop

at its distal end and the endosteal guide 45 may include some sort of hook at its distal end, or the

end of the endosteal guide 45 may include some sort of loop at its distal end and the end of the

wire suture 55 may include some sort of hook at its distal end, etc. Alternatively, endosteal

guide 45 may be cannulated so as to permit grippers or a hook to be passed down the interior of

the endosteal guide 45 to grab the distal end of the wire suture 55. In any case, the endosteal

guide 45 and wire suture 55 are configured so as to permit the distal end of the endosteal guide

45 to connect up with the distal end of the wire suture 55.

[10043][0034] Referring next to Fig. 9, once the distal end of endosteal guide 45 is connected

with the distal end of wire suture 55, a cannulated screw 60 is slid down the wire suture 55 and

screwed into angled tunnel 50. This can be done by using a cannulated driver (not shown in Fig.

9) of the sort well known in the art. Cannulated screw 60 is screwed into femur 10 until the head

of the screw 60 contacts the periosteum so that the head of the screw 60 does not stand proud

above the bone. Cannulated screw 60 essentially acts as a sort of liner for angled tunnel 50,

providing a bearing surface for wire suture 55 so as to prevent the wire suture 55 from cutting

into femur 15 as the wire suture 55 is moved through angled tunnel 50, as will hereinafter be

discussed.

100441 [0035] Next, and looking now at Fig. 10, drill guide 40 is withdrawn from the surgical

site. As this occurs, the drill guide's endosteal guide 45 of drill guide 40 is removed from bone

Page 7

tunnel 20, causing wire suture 55 to be pulled down femoral tunnel 35, across the interior of the knee joint, down tibial tunnel 30 and then out the bottom end of the tibial tunnel 30.

<u>[0045][0036]</u> Then, as shown in Fig. 11, the graft ligament 25, e.g., one or more strands of a hamstring tendon, is attached to the distal end of wire suture 55. In the case where the graft ligament <u>25</u> comprises one or more strands of hamstring tendon, such attachment can be easily accomplished by forming a loop at the lower end of wire suture 55 and then passing the tendon strand(s) through the loop, in the manner shown in Fig. 11.

[0037] Next, as shown in Fig. 12, wire suture 55 is pulled proximally so as to tow graft ligament 25 up bone tunnel 20 until the distal end of the graft ligament 25 is positioned in femoral tunnel 35 and the proximal end of the graft ligament 25 is positioned in tibial tunnel 30. Then a crimp 65 is attached to wire suture 55 adjacent to cannulated screw 60 so that engagement of crimp 65 with cannulated screw 60 will prevent the wire suture 55 from being pulled back in the tibial direction. Next, graft ligament 25 is tensioned by pulling the proximal end of the graft ligament 25 in the proximal direction, and then the. The proximal end of the graft ligament 25 is then fixed to the tibia in ways well known in the art.

<u>[4047][0038]</u> If desired, the proximal end of cannulated screw 60 can have a modest recess therein so as to receive crimp 65 when the ligament is tensioned, whereby to keep the crimp from standing proud above the femur. Alternatively, the bone screw <u>60</u> can be eliminated, utilizing only the crimp <u>65</u> on the wire <u>suture 55</u> to prevent distal migration. Finally, the end of wire suture <u>55</u> proximal to crimp <u>65</u> is trimmed off, thus effectively completing the ACL reconstruction procedure.

[0048][0039] In the foregoing description, and in Figs. 11 and 12, graft ligament 25 was discussed in the context of comprising one or more strands of hamstring tendon. However, it should also be appreciated that graft ligament 25 may comprise other constructs as well. Thus,

for example, graft ligament 25 may comprise one or more bone blocks attached to a tendon (e.g., a so-called "bone-tendon-bone" graft), or a totally artificial prosthesis, etc.

<u>[4049][0040]</u> By way of example but not limitation, Figs. 13 and 14 shown show a bone-tendon-bone graft being towed into position within bone tunnel 20. In this case, wire suture 55 may be attached to the leading bone block 70, e.g., by passing the wire suture <u>55</u> through one or more holes formed in the bone block <u>70</u> and making it fast. If desired, a second wire suture <u>75</u> may be attached to the trailing bone block <u>80</u> so as to fasten the trailing end of the ligament graft to tibia 10 in ways well known in the art.

[0050][0041] It should be appreciated that the aforementioned procedure and apparatus may be modified without departing from the scope of the present invention.

[0042] Thus, for example, and looking now at Fig. 15, in one alternative reconstruction procedure, after angled tunnel 50 has been drilled, drill guide 40 is completely removed from the surgical site. Then a wire suture 55, preferably having a loop 83 at its distal end, is passed down angled tunnel 50 until the wire suture 55 begins to pass into the top end of femoral tunnel 35.

H00521[0043] Next, a suture grasper 85 passed up bone tunnel 20 until the distal end of the suture grasper 85 is positioned adjacent to the distal end of wire suture 55, whereupon the. The suture grasper ean be-85 is then used to pick up the distal end of wire suture 55 and draw the wire suture 55 down femoral tunnel 35, across the interior of the knee joint, down tibial tunnel 30, and finally out the front of tibia 10, as shown in Fig. 16. The cannulated screw 60 can then be mounted on the proximal end of wire suture 55, the cannulated screw 60 screwed down into the periosteum, the distal end of wire suture 55 fastened to the graft ligament 25, the wire suture 55 retracted so as to tow the graft ligament 25 up bone tunnel 20 and into position, the crimp 65 applied to the wire suture 55, and then the graft ligament 25 made fast to the tibia, all in substantially the same manner as previously described.

[0053][0044] And in In another alternative reconstruction procedure, cannulated screw 60 and crimp 65 may be replaced by an alternative construction. More particularly, and looking now at Figs. 17-19, a cannulated screw 90 and an associated threaded collet 95 are shown. Cannulated screw 90 has the usual external threads 100 (Fig. 17) and central bore 105 (Fig. 18); Cannulated screw 90, however, cannulated screw 90 also has a threaded counterbore 110 (Fig. 18) terminating, intermediate the cannulated screw 90, in a smooth tapered section 115. Correspondingly, threaded collet 95 comprises a leading tapered section 120 (Fig. 17), a trailing threaded section 125, and a central bore 130 (Fig. 18). The collet's leading tapered section 120 of collet 95 is slit, whereby the threaded tapered section 120 can be forced radially inward so as to close down the diameter of the eollet's-central bore 130.—Cannulated of collet 95. Depicted in Fig. 20, cannulated screw 90 and threaded collet 95 are sized so that the threaded collet 95 may be lightly screwed into the rear of cannulated screw 90 and wire suture 55 threaded therethrough: thereafter. Thereafter, threaded collet 95 may be screwed further into cannulated screw 90 so that the nose of the threaded collet 95 will close down on wire suture 55, whereby to elampthereby clamping the wire suture 55 to the threaded colletcollet 95 and, hence, to the cannulated screw 90.

[0054][0045] Figs. 17 and 21 show drivers for advancing and, alternatively, retracting cannulated screw 90 and threaded collet 95. More particularly, a first cannulated driver 135 is provided for turning cannulated screw 90. A second cannulated driver 140 is provided for turning threaded collet 95. Second cannulated driver 140 may be positioned within first cannulated driver 135, as will hereinafter be discussed.

[0046] In use, after wire suture 55 had been passed through angled tunnel 50, cannulated screw 90, threaded collet 95, second driver 140 and first driver 135 are mounted on the proximal end of wire suture 55 (Fig. 17). Next, the components are manipulated so that threaded collet 95

is lightly threaded into cannulated screw 90, second driver 140 is engaged with threaded collet 95, and first driver 135 is engaged with cannulated screw 90. Then the aforementioned assembly is moved down wire suture 55 until the leading tip of cannulated screw 90 engages the femur.

whereupon first driver 135 is used to drive the cannulated screw into the femur.

[0056][0047] At this point, cannulated screw 90 will line lines angled tunnel 50, but wire suture

55 will-be-is free to move relative to cannulated screw 90 and, hence, the patient's anatomy.

Wire suture 55 is then used in the manner previously described to pick up graft ligament 25 and

tow the graft ligament 25 back up into position within the bone tunnel. Once the graft ligament

25 is in position, wire suture 55 is made fast by screwing threaded collet 95 further into

cannulated screw 90, using second driver 140, until the leading tip of the threaded collet 95

closes down on wire suture 55, whereupon the wire suture will be 55 is clamped to the threaded

collet 95 and, hence, the cannulated screw 90. Once this has been achieved, drivers 135 and 140

may be removed, the excess wire 55 trimmed away, and the tibial side of the graft secured.

[0057][0048] It should also be appreciated that the procedure and apparatus described above

may be used for purposes other than an ACL repair, e.g., they may be used to repair other

ligaments, the apparatus may be used in other types of surgical procedures such as trauma, spine,

etc.

10058100491 Also, the wire suture may be braided polyethylene or monofilament suture; and

the cannulated screw and/or threaded collet may be plastic or even reabsorbable.

Page 11

Abstract

Apparatus is disclosed for reconstructing a ligament, the apparatus including a bone

fixation element having a central bore and adapted for positioning in a bone tunnel, a flexible

filament configured for holding a graft ligament, and a crimp configured for attachment to the

flexible filament so as to prevent movement of the graft ligament in the bone tunnel.

W:\13447\46\DFW0000014096V001.doc

Page 12

Serial No. 09/970,559
Amendment Dated May 31, 2005
Reply to Office Action dated December 1, 2004
Annotated Sheet

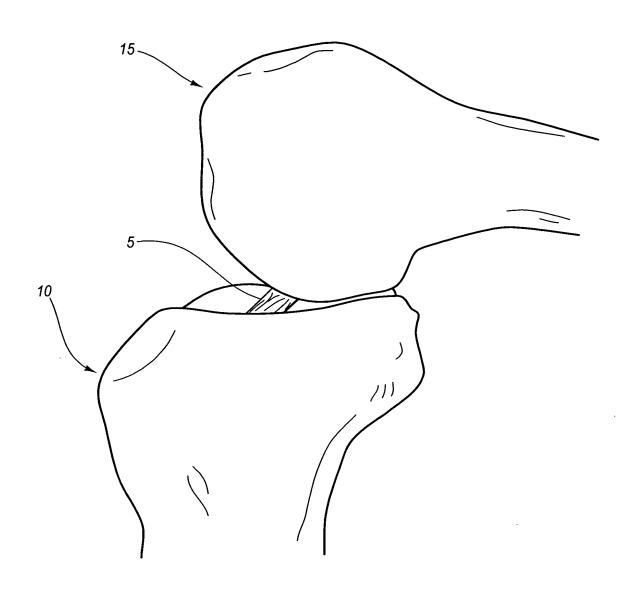


Fig. 1

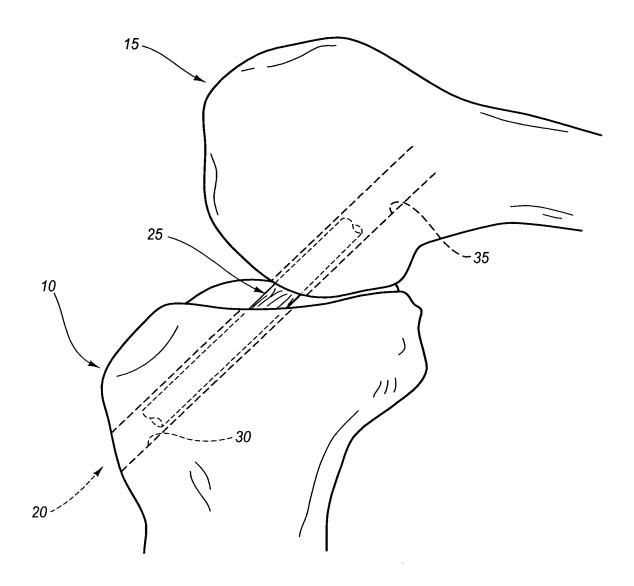


Fig. 2

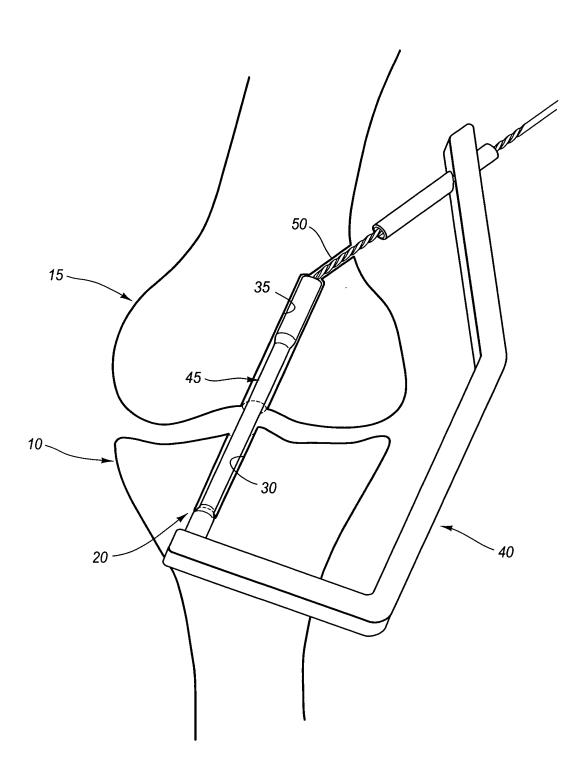


Fig. 3

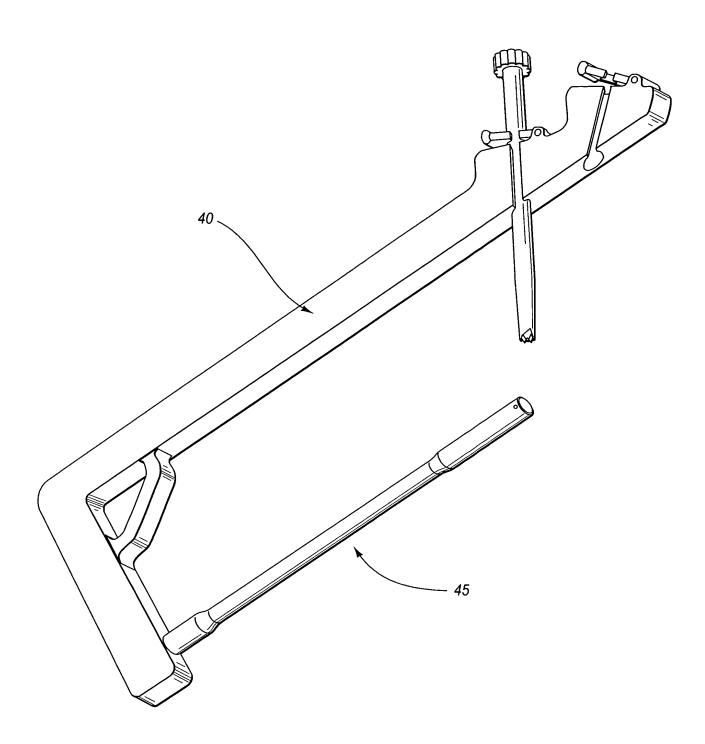
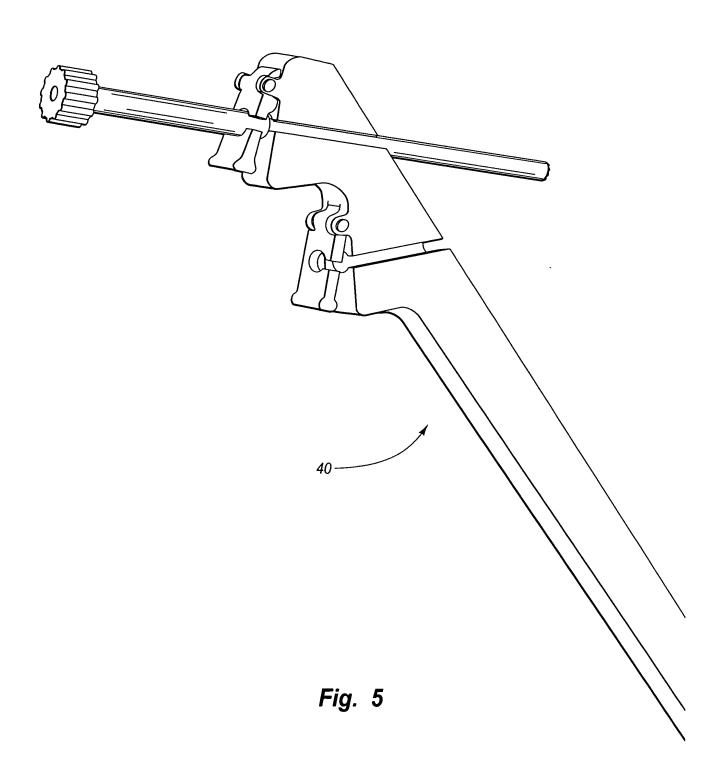


Fig. 4



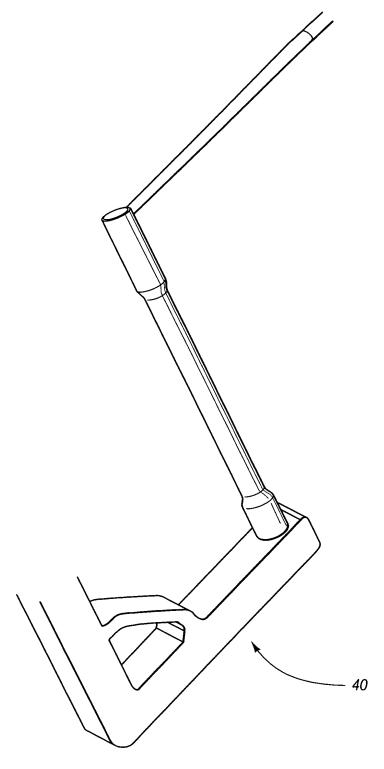


Fig. 6

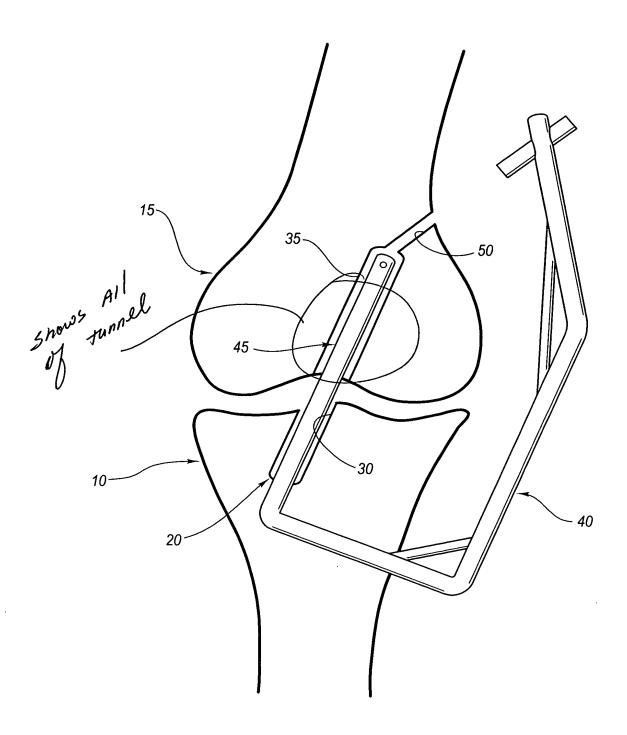


Fig. 7

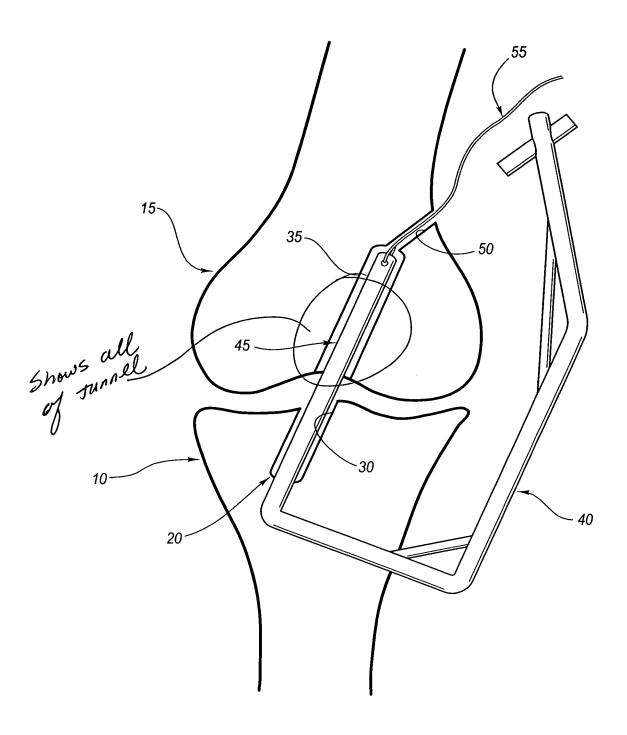


Fig. 8



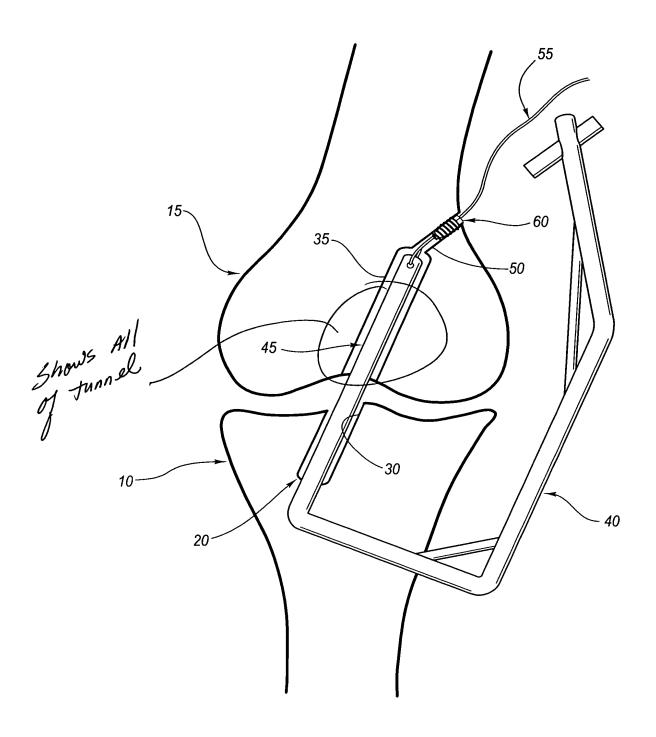


Fig. 9

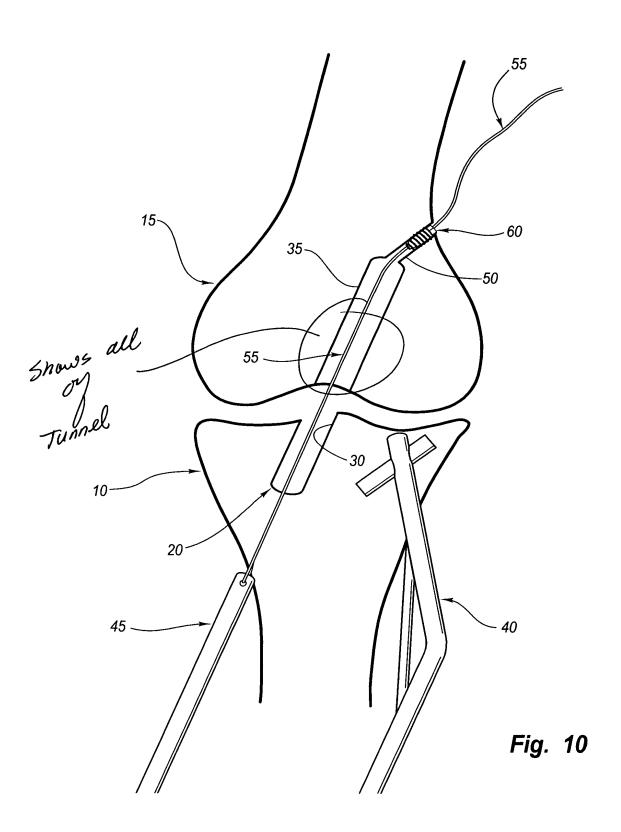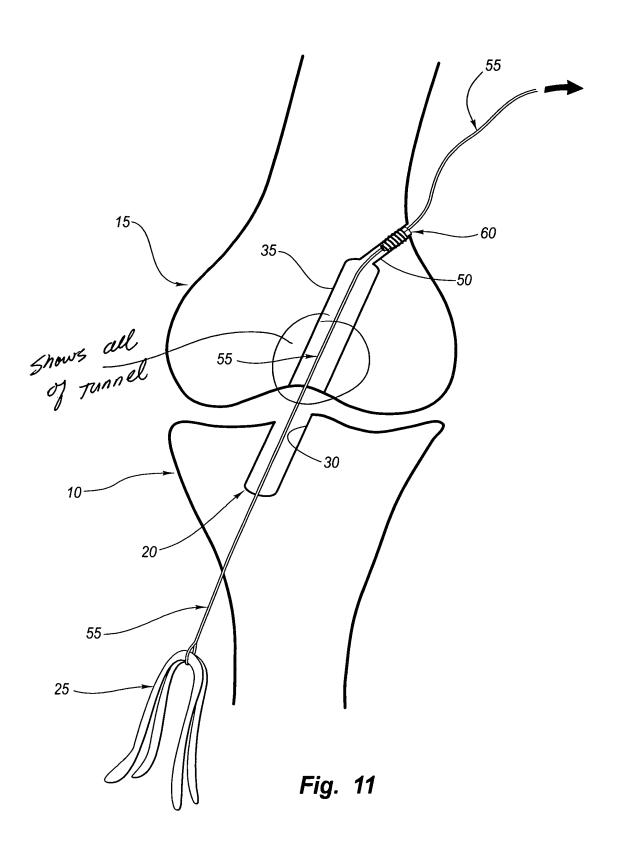

11 / 20



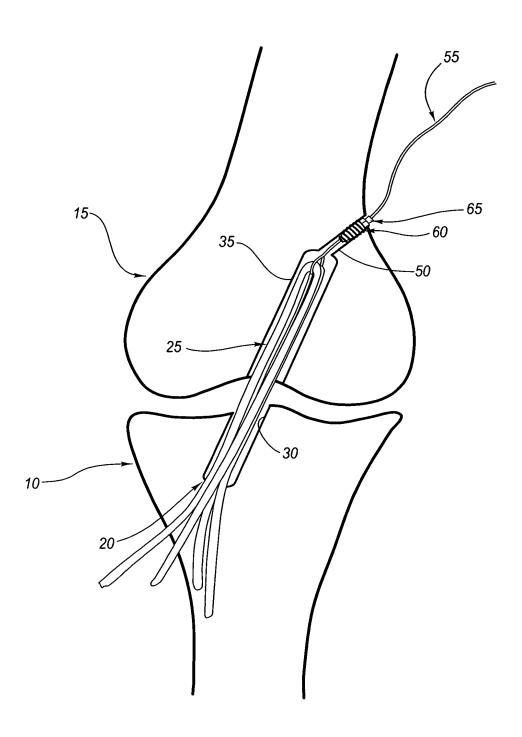
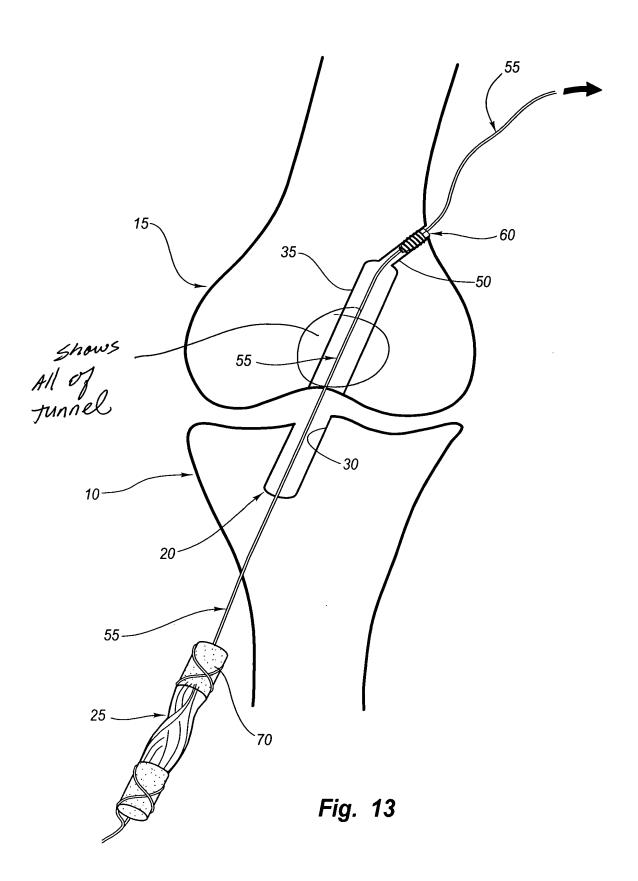


Fig. 12



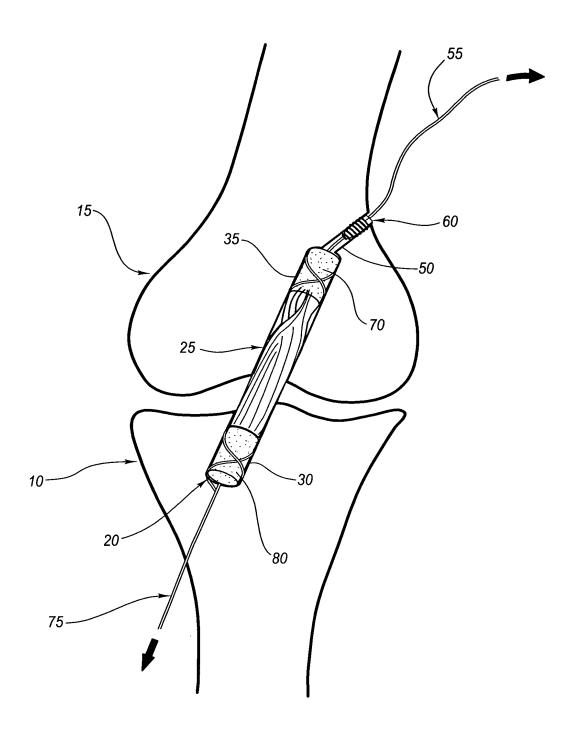


Fig. 14

15/20

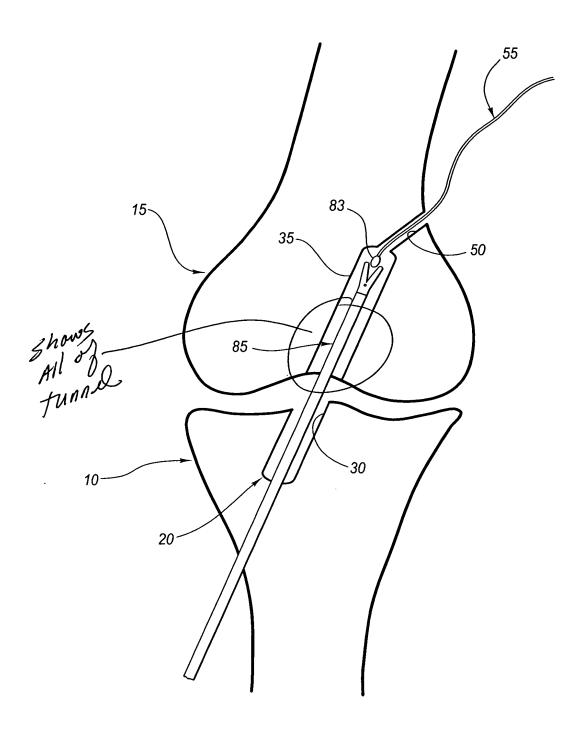
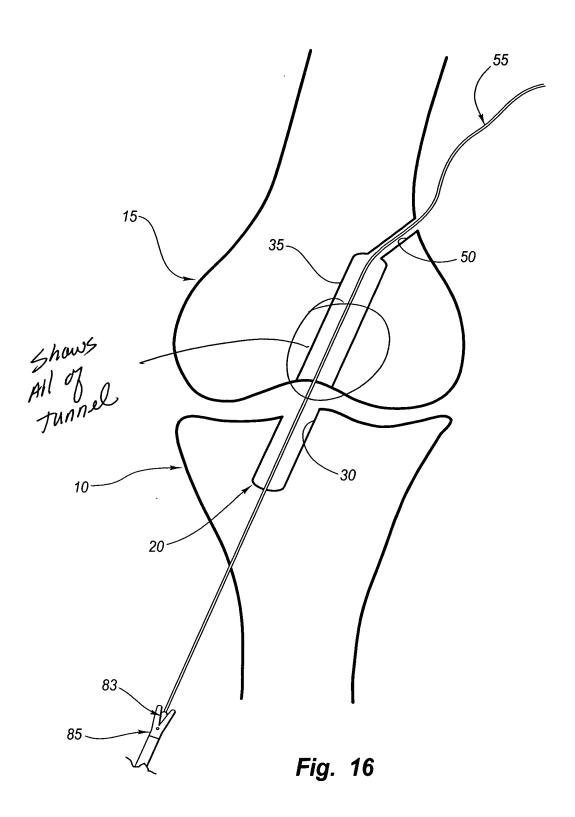
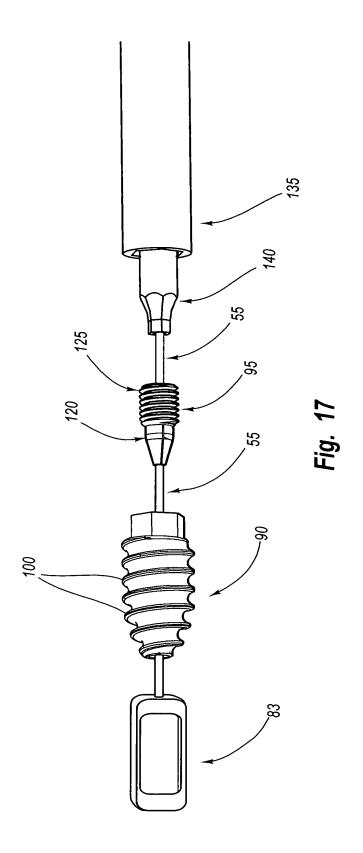


Fig. 15

16/20





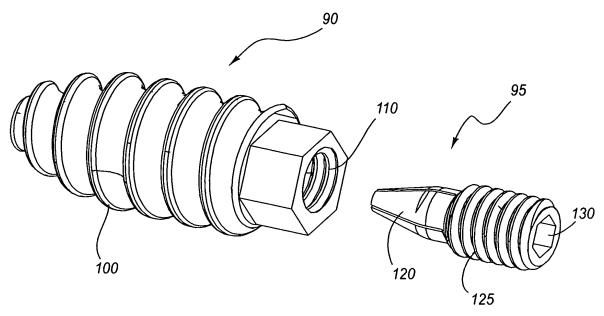


Fig. 18

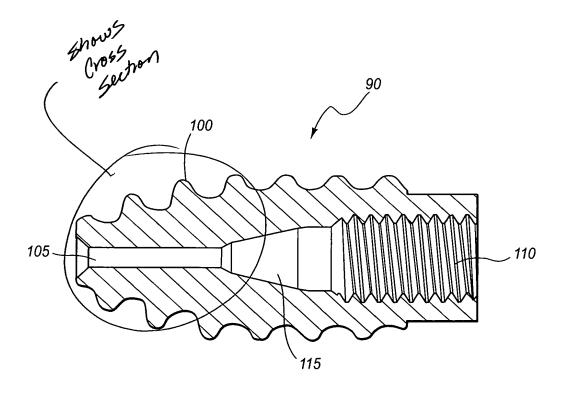


Fig. 19

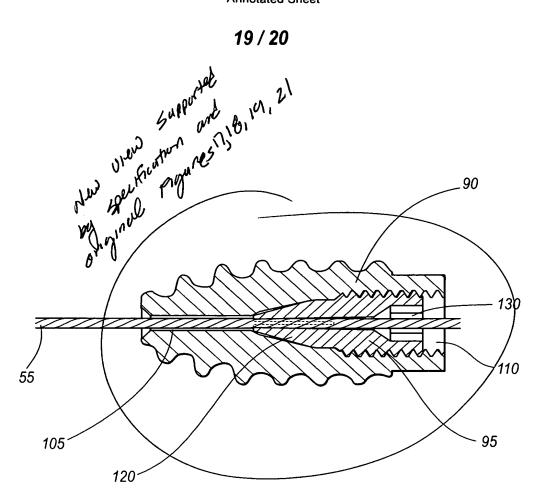


Fig. 20

Serial No. 09/970,559
Amendment Dated May 31, 2005
Reply to Office Action dated December 1, 2004
Annotated Sheet

